

Ready-2-Use



The 2 in 1 system with disposable cleanroom mop and integrated disinfectant

Hydroflex

PurMop® Ready-2-Use

The innovative 2 in 1 solution: PurMop® disposable mops with integrated disinfectant from InSpec™. Select your individual solution!

The completely new product concept from Hydroflex and Redditch Medical uses a special burst pouch technology to provide a disposable mop with a full range of cleaning and disinfecting agents specially designed for cleanrooms.

The disposable mop and disinfectant are delivered in one packaging unit but are separated through a burstable compartment. For use simply burst the pouch, distribute the fluid and after a short waiting time the disinfectant will evenly and sufficiently saturate the mops. Thanks to the separation of the two components a long shelf life for mop and fluid can be ensured without reducing the effectiveness of the disinfectant over time.

Your advantages:



Safe

Safe means that you always get a new and fully functional mop which will be soaked with the correct amount of disinfectant. There is no need to validate an additional trolley system.



Simple


Simple means that you only have to combine the fluid with the mop and after a short waiting time you can start with cleaning. Neither cleaning trolleys nor the storage of additional disinfectants is necessary.



Innovative

Rethinking disinfection of critical cleanroom areas! With the PurMop® InSpec™ Ready-2-Use concept we are offering a completely unique and innovative solution for efficient cleanroom disinfection. As a 2 in 1 On-Demand packaging, this solution is simple, fast, safe and eliminates the need for additional storage of disinfectants or a trolley system.

IDEAL FOR SMALL
STERILE CLEANROOMS! ■

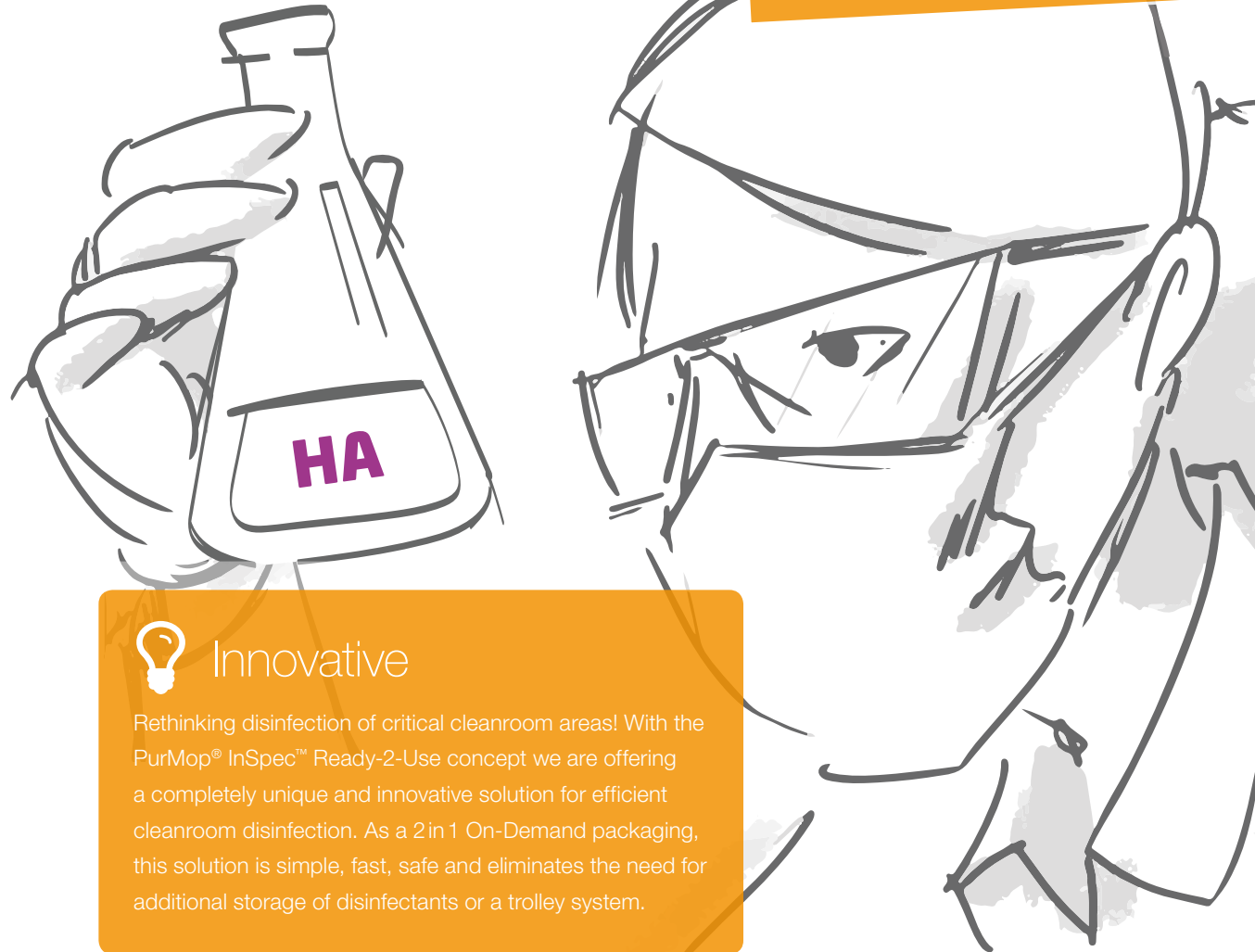

1

Single-Use Mop



EF40:
Polyester Microfiber Mop
for cleaning and disinfection
of floors, walls and ceilings



EFB40:
Polyester Microfiber Mop with
increased cleaning ability for
floors with higher contamination



EC20:
Polyester Knitted Mop
for disinfection of isolators
and workbenches

Mix & Clean



2 in 1
system

Disinfectant



QT
Quaternary ammonium
compound



AN
Diamine



HA
Hypochlorous acid



OX
Peracetic acid
and hydrogen peroxide



N10
Alkaline detergent

Product	Bactericide	Fungicide	Sporicide
QT	✓	✓	-
AN	✓	✓	-
HA	✓	✓	✓
OX	✓	✓	✓

2

How it works:




Activate: Mix the fluid with the mop
by bursting the membrane and shake it
several times.




2 MINUTES

Saturate: Wait 2 minutes until the mop
is evenly saturated

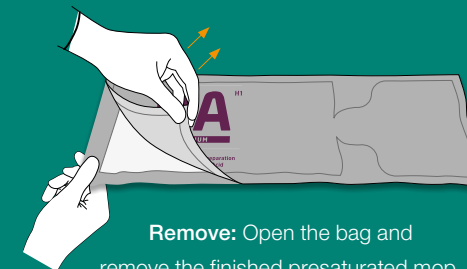

Remove: Open the bag and
remove the finished presaturated mop.

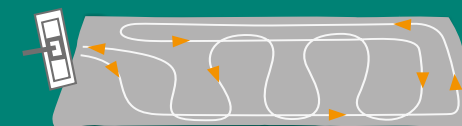

Wipe: Start wiping. **Finished!**



PurMop® Safe-Use

Single-Use Mop



InSpec™

Detergent and disinfectant



HYDROFLEX OHG Cleanroom Hygiene

Am Weidenhäuser Bahnhof 8 ▶ 35075 Gladenbach
Germany ▶ T +49 6462 91598-0 ▶ F +49 6462 91598-20
hydroflex-group.com



REDDITCH MEDICAL ▶ Discovery 2

2 William Armstrong Way ▶ NETPark ▶ Co Durham
Sedgefield ▶ TS21 3FD ▶ UK ▶ T +44 1527 830942
redditchmedical.com

Contact Hydroflex for further information or local distribution
info@hydroflex-group.com